




Laryngoscopes

Fiber Optic Reusable Laryngoscopes


Macintosh

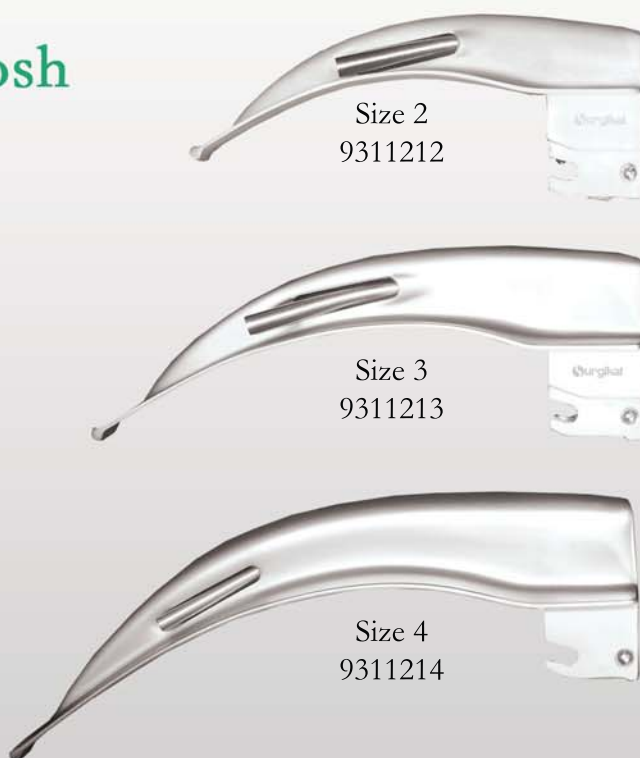

What should I know about Laryngoscope blade shapes?

We offer the two standard blade shapes. Macintosh is the most widely used of the curved blade type and is the preferred shape for most adult intubations. Miller is straight, and designed to be longer and narrower than other blades. Most pediatric and neo-natal intubations are performed with Miller.

Miller

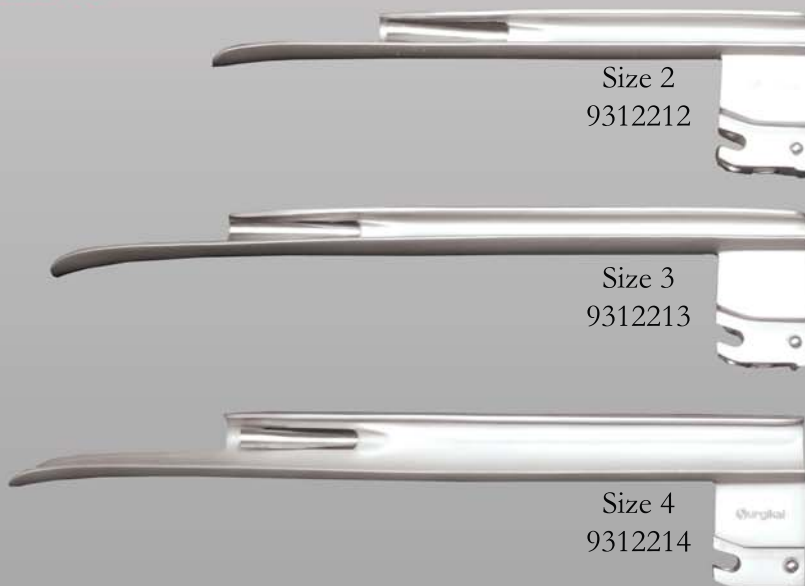

Fiber Optic Disposable Laryngoscopes

Macintosh



What should I know about Fiber Optic lighting systems?


Fiber Optic transmits the illumination from our xenon lamp in the handle through the fiber optic tube in the blade. The fiber optic tube is an excellent conduit of the light and this system is used by many anesthesiologists. Fiber Optic systems, handles and blades, are indicated by green marking.

Miller



Conventional Reusable Laryngoscopes

Macintosh

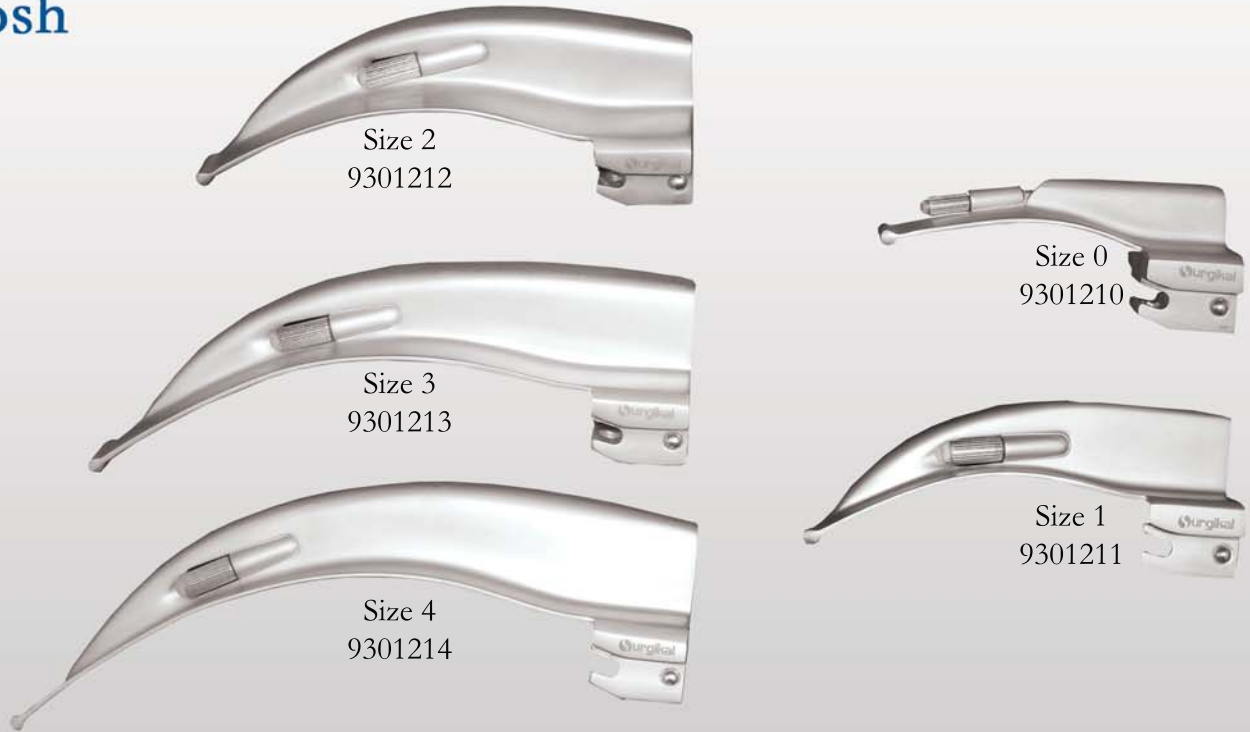

What should I know about Conventional lighting system?


Conventional lighting has the lamp in the blade. The light is activated when the contact point in the blade meets the contact point in the handle. PMD uses a xenon lamp in its blades to provide a strong white light. The conventional system is growing rapidly in popularity.

Miller



Conventional Disposable Laryngoscopes

Macintosh

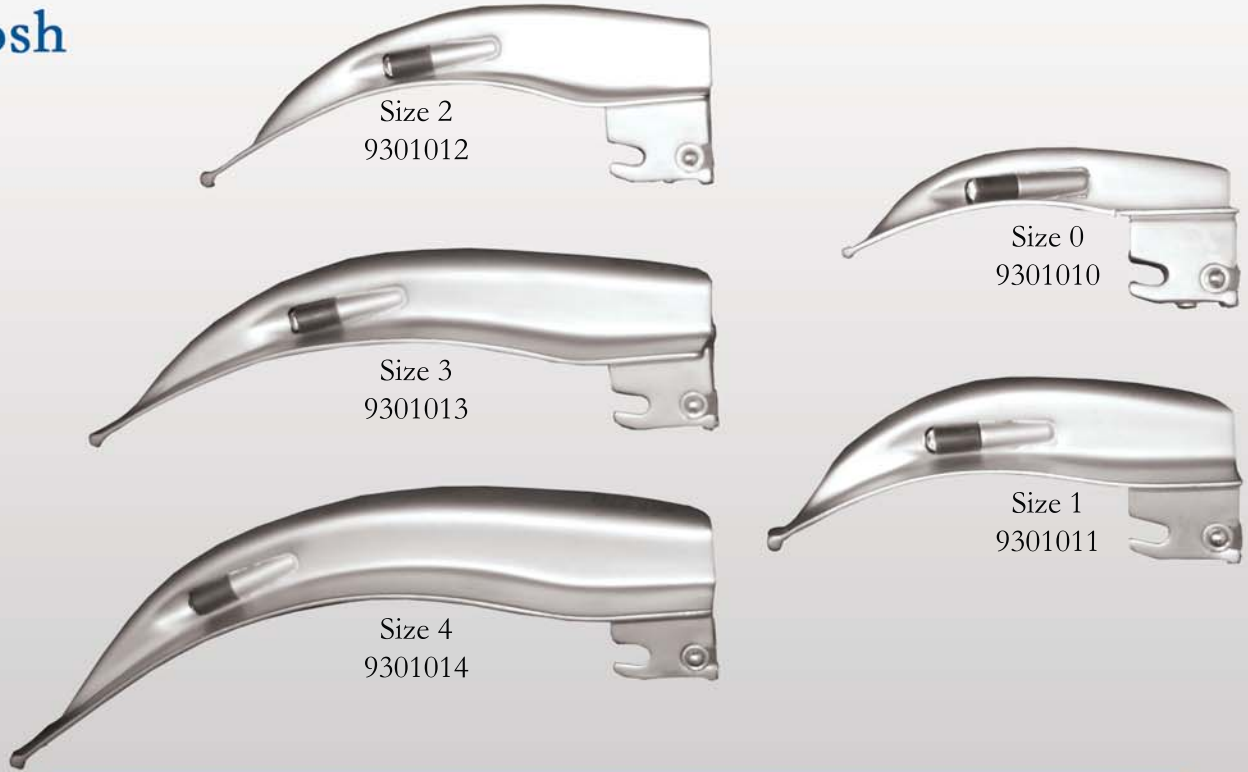

What should I know about Disposable Laryngoscopes blades?

Disposable or single-use blades have become the blade of choice in many institutions. They provide cost advantages, save re-processing time and expense, and help prevent cross contamination. PMD has one of the broadest range of disposable surgical grade stainless steel blades available.

Miller



Specialty Laryngoscope

Laryngoscope Handles



Medi-Scope

High quality stainless steel ENT single patient use system.

Strong conventional lighting system.

9390203 - Mac 3

9390204 - Mac 4



Conventional Handles

Fiber Optic Handles

Large	Medium	Small	Stubby
9310003	9310002	9310001	9310000
9300003	9300002	9300001	9300000

What should I know about PMD?

Precision Medical Devices (PMD) is a leading manufacturer of high quality, great value laryngoscope blades and handles with a wide selection of the two most used blade shapes: Macintosh and Miller. Special sizes such as the Miller 00 are available. Our handles are manufactured of nickel coated brass with balanced heft and sure grip surface.

Our products use a xenon lighting system providing strong, white light. Our Fiber Optic system is an excellent conduit of the illumination and this system is a traditional favorite of many anesthesiologists. For those who prefer Conventional lighting, our fastest growing category, we offer a complete range of shapes and sizes.

Disposable or single-use blades have become the blade of choice for many institutions. They provide cost advantages, save re-processing time and expense, and help prevent cross contamination, a key objective of Joint Commission. We manufacture our disposable blades only from surgical grade stainless steel.

All Precision Medical laryngoscope blades and handles are manufactured under strict supervision to our exacting quality standards.

Distributed by:

(732) 447-2587



www.PMDmfg.com